

OBSTETRIC PRODUCTS


RELIABLE AND COST EFFECTIVE



Amniotic Membrane Perforator

The Amniotic Membrane Perforator provides a safe and comfortable method for rupturing the amniotic membrane. The Amniotic Membrane Perforator is thoroughly inspected and rigorous quality control systems ensure a high level of quality.

- Rounded and protected tips to perforate the amniotic membrane without endangering the mother or the baby.
- Disposable, sterilized, and individually packed for hygienic application.
- Slipping is prevented thanks to the ergonomically shaped sides and non-slip surface.
- The amniotic membrane perforator is 26.5 cm long.

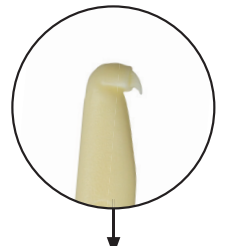

Umbilical Cord Clamp

This precision Umbilical Cord Clamp is hygienic, safe and reliable to use. After the birth, the umbilical cord can be clamped to prevent blood loss in the newborn baby.

- Safe disposable product, easy to use.
- The non-slip ends of the umbilical cord clamp prevent slipping during use.
- Seals securely and exerts even pressure regardless of the diameter and contraction of the umbilical cord.
- The guard on the hinge holds the umbilical cord between the toothed surfaces and fixes this permanently.



Umbilical Cord Clamp Clipper

The Umbilical Cord Clamp Clipper was developed for severing the specially constructed hinge of the Umbilical Cord Clamp and removing this from the baby.

- Safe and simple handling.
- The concealed cutting surface protects the baby from injury.
- Disposable usage prevents the risk of cross-contamination.



OBSTETRIC PRODUCTS

RELIABLE AND COST EFFECTIVE



Urine Specimen Bag

The PVC-free Urine Specimen Bags are the first choice when it comes to convenience and protection for delicate skin. The sophisticated system for urine collection can be used for babies and children of both sexes.

- Different shapes and sizes to handle almost all anatomical differences.
- Standard or perforated adhesive surface (both latex-free).
- The double-chamber construction with a non return-flow valve keeps the urine away from the skin and prevents skin irritation.
- Adhesive surface can be folded over and stuck down to send the urine sample to the laboratory.
- The horizontally split protective film on the adhesive surface enables easy application.

200 ml Pediatric




100 ml Newborn



10 ml Preemie



Porous adhesive

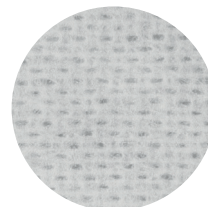

Standard adhesive



You can find more information at:
www.ardomedical.com/obstetrics

99.04.341 Ardo OB Flyer en V01 10.2019.pdf